# WelchAlllyn®

## English

### 719 Series Lithium Ion Handle

Material No. 723771 Rev. B  
© 2015 Welch Alllyn, Inc. Revision date: 2015-08

Regulatory Affairs Representative  
Welch Alllyn Limited  
Navan Business Park  
Dulinn Road  
Navan, County Meath,  
Republic of Ireland

### INTRODUCTION

Thank you for purchasing the Welch Alllyn 719 Series Lithium Ion Power Handle. The handle is intended to be used as a power source for all Welch Alllyn standard instrument heads. The 719 series represents a new generation of handles from Welch Alllyn. Taking advantage of new battery technology, the 719 series provides you double the on-time in a smaller and lighter weight package. Smarter battery technology means that "memory" is no longer an issue. What's more, the 719 reminds you to recharge with its low battery indicator. The soft rheostat section provides a sure grip and the overall design is ergonomically sound. The Charging Base and Charger charges the handle in 14 hours.

**Before using this product, please read the following warnings and instructions.**

### SAFETY AND WARNINGS

 Attention: Read Operating Manual for Cautions and instructions for use.

- POSSIBLE EXPLOSION HAZARD. Equipment not suitable for use in the presence of a flammable anesthetic.
- Use the 719 Series Handle only with Welch Alllyn approved products and chargers.
- Welch Alllyn p/n 71960 (Battery Section) is only suitable for use with the 719 series handle.
- Do not disassemble or modify the battery pack. There are no serviceable parts inside.
- Do not attempt to directly solder to the battery pack.
- Do not attempt to connect the positive and negative terminals to each other, nor to any other device.
- Do not charge battery pack in any environment where the temperature may exceed 40°C (104° F) or fall below 0°C (32° F).
- Do not dispose of the battery section in fire or water. Do not puncture, incinerate or subject to strong impact or shock.
- Do not immerse the battery pack in any liquid.
- Do not strike or throw the battery pack.
- Do not use a battery pack which appears to be deformed, leaking, corroded, or is otherwise irregular.
- Handle a damaged or leaking battery with extreme care. If you come into contact with the electrolyte, wash the exposed area with soap and water. If it contacts the eye, flush the eye with water for 15 minutes. Seek medical attention immediately.
- When discarding a battery pack, contact your local waste disposal service provider regarding local restrictions on the disposal or recycling of batteries.
- To avoid the risk of fire, burn, injury, or damage to your battery pack or lamp, do not allow a metal object to touch across the battery contacts.
- This product complies with current required standards for electromagnetic interference and should not present problems to other equipment or be affected by other devices. As a precaution, avoid using this device in close proximity to other equipment.

FOR BEST PERFORMANCE:

- Do not expose battery pack to storage temperatures above 50°C (122°F) or below -20°C (-4°F).
- Do not discharge battery pack if the ambient temperature is above 60°C (140°F) or below -20°C (-4°F).

### Parts List

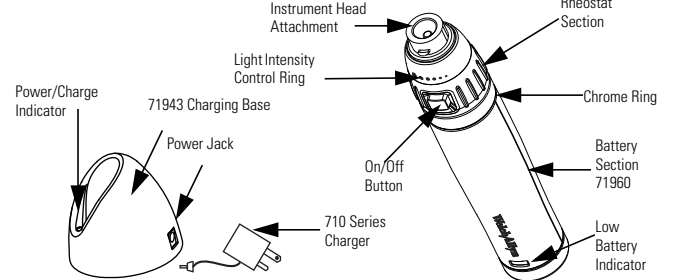
Figure 1

### OPERATING INSTRUCTIONS

- Connect instrument head to the handle.
- Depress On/Off Button on Rheostat Section and rotate the Light Intensity Control Ring clockwise (CW). To intensify light, continue to rotate until stop is reached (Figure 1).
- To turn off, rotate counter clockwise (CCW). Button will click when fully off. Turn fully off after each use to ensure maximum on-time from the battery.

### CHARGING INSTRUCTIONS

- Plug the Charger into a wall outlet and the Power Jack on the Charging Base. The Power/Charge Indicator will illuminate (green) when power is applied (Figure 1).
- Handle the handle the Charging Base. The Power/Charge Indicator will turn to amber during charge. Charge for 14 continuous hours to charge a dead battery (Figure 2).
- When charging is complete, the Power/Charge Indicator will return to green. The handle may remain in the Charging Base if desired.
- The Battery may be charged after partial discharge without negatively affecting overall battery memory (unlike Ni-Cd cells).
- Diagnostic Instruments cannot be used while charging.

**CAUTION** The Battery Section is not shipped fully charged. For maximum performance charge the Battery Section for 14 hours when new or after a long period of inactivity.

NOTE: The handle is only compatible with the 71943 charging base, 71950 AC charger module, 7114X and 7111X (with adapter) desk chargers.


### LOW BATTERY INDICATOR

- The Low Battery Indicator will illuminate several minutes before handle shuts off depending on the age of the battery (Figure 1).
- Once the battery is fully discharged, the Low Battery Indicator may flicker.

### BATTERY REPLACEMENT

#### Removing the Battery Section:

- Hold the instrument handle firmly by the Battery Section
- Turn the Rheostat Section counterclockwise until it detaches from the Battery Section (Figure 3).
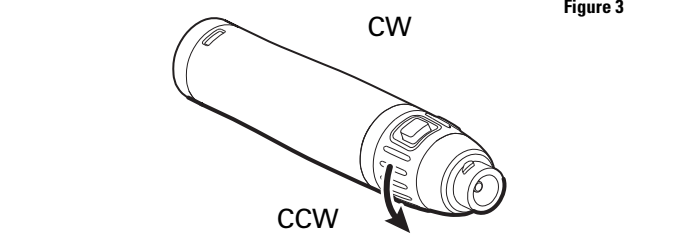
#### Replacing a New Battery Section:

- Insert the Battery Section into the Battery Section
- To connect the two sections, turn the Rheostat Section clockwise until tight.
- Verify the Battery Section is secure on the Battery Section of the instrument handle.

### CLEANING AND DISINFECTION

Wipe the handle and charger with an appropriate health care low- or intermediate-level cleaner/disinfecting wipe that incorporates either a 1:10 sodium hypochlorite (bleach) solution or isopropyl alcohol as the active disinfection ingredient. Follow wipe manufacturer's instructions for appropriate use, contact times and applicable warnings and precautions.

Do not excessively soak/saturate the handle assembly or charger.

Do not immerse handle assembly or charger in any solution.

Do not sterilize the 719 Series power handle, charging base, or charger.

### Battery Recycling

Welch Alllyn now employs the services of an agency which can disassemble and recycle all components of rechargeable batteries so that nothing gets landfilled or incinerated. Should you wish to dispose of your expended Welch Alllyn rechargeable battery via recycling (in the U.S. and Canada only), please send to:

Welch Alllyn Canada  
RECYCLE BATTERY  
4341 State Street Road  
Mississauga, Ontario  
M5S 1B3  
CANADA (42 1V4)

For locations outside the U.S. and Canada consult your local recycling authority for correct disposal of batteries.

### SPECIFICATIONS

(Specifications and design may change without notice)

#### Handle

**Model Number:** 71911  
**Dimensions:** 5.4" L x 1.125" Dia. (Max.) (137mm x 28.6 mm)  
**Weight:** 4.5 oz. (128 g)  
**Typical Charge Time:** 14 hours (fully discharged)  
**Typical Discharge Time:** 110 minutes (for new battery)  
**Conforms with:** IEC/UL/CSA/EN 60601-1 when used with Charging Base 71943 and 710 Series Charger

#### Charging Base

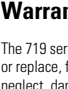
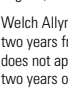
**Model Number:** 71943  
**Conforms with:** IEC/UL/CSA/EN 60601-1 when used with 71911 Handle & 710 Series Charger

#### Charger

**Model Number** Input  
71030 100V ~ 50/60Hz, 50mA  
71032 230V ~ 50/60Hz, 25mA  
71034 230V ~ 50/60Hz, 25mA  
71036 240V ~ 50Hz, 25mA  
**Conforms with:** All Charger Model numbers conform with IEC/UL/CSA/EN 60601-1

#### Continuous Operation

#### SYMBOLS

	<b>Environment (All System Components)</b> Temperature limitation Operating: 59°F (10°C)-104°F (49°C) Transport/storage: -4°F (-20°C) - 120°F (55°C)
	Humidity limitation Operating: 30% - 90% Transport/storage: 10% - 95%
	Atmospheric pressure limitation 500 hPa - 1060 hPa
	Attention. Read Operating Manual for Cautions and Instructions for Use.
	Type BF applied part
	Class II equipment
	Reorder number
	Product Identifier
	Manufacturer & Date of Manufacture (YYYY-MM-DD)
	Professional use only
	Authorized Representative in the European Community
	Global Trade Item Number
	Consult operating instructions
	Batch code
	Equipment not protected against the ingress of water.

### SERVICE AND WARRANTY

#### Service

If the product fails to function properly or if assistance, service, or spare parts are required, contact the nearest Welch Alllyn Technical Support Center.

**Welch Alllyn Technical Support:**  
http://www.welchalllyn.com/about/company/locations.htm

### Warranty

The 719 series Power Handle is guaranteed by Welch Alllyn against all manufacturing defects. Welch Alllyn will repair or replace, free of charge, any parts of its own manufacture proven to be defective through causes other than misuse, neglect, damage in shipment, or normal wear.

Welch Alllyn warrants that the Welch Alllyn 71960 Battery Section will perform to original specifications during the two years from the date of manufacture, when used with the Welch Alllyn 719 series Power Handle. This warranty does not apply to the use of the 71960 Battery in other products. A defective battery will be replaced if it fails within two years of the date of manufacture.

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